

Practical Batch Record- Wet Granulation- V04

Page 1 of 19

DoE Trial No: CP (Low Shear)

Practical group No: 11

Batch Record Summary

1	PRE REQUISITES	2
2	PROTOCOL INFORMATION	3
3	DISPENSING	4
4	STEP 1: DRY MIXING	6
5	STEP 2: GRANULATION	7
6	STEP 3: DRYING	8
7	STEP 4: SCREENING	9
8	STEP 5: LUBRICATION	9
9	STEP 6: TABLETTING	12
10	STEP 7: TABLETTING GAMLEN.....	15
11	GENERAL COMMENTS ON THE BATCH.....	16
12	CLEAN LABELS.....	17

Signature prior to the trials Date/sign	Project Manager	QA Manager	Scientists	Signature during the trials Date/sign	Scientist1	Scientist2	Check signature Date/sign	QA Manager
VN 07 APR 2022.	VS 07 APR 2022.	Yushu		Turk	KS		Lilly 07 APR 2022.	

7/4/22

7/4/22

7/9/22

Practical Batch Record- Wet Granulation- V04

DoE Trial No:

CP1

Practical group No:

II

1 GMP PRE REQUISITES

1.1 Signatures

Before starting, read the whole document and populate the header information on every page and proceed to the signatures at the right locations for the beginning of the trial. Each page will also have to be signed during the trial and after the trial by the required persons.

1.2 Raw data recording

THE BATCH RECORD HAS TO BE FILLED WITH BLUE INK PENS AND THE QA VERIFICATION HAS TO BE DONE WITH A RED PEN

If you make a mistake when you record the raw data, you have to be clearly cross it (the mistake needs to be still legible), you put the date and you sign (if you have room you have even to explain the mistake)

Start date and time must all be completed.

1.3 Clean labels

You must start the trial or testing only if the manufacturing or testing equipment are CLEAN and have a "Clean" label partially stuck onto one part of the equipment. That "Clean" label must contain the following information:

- Equipment name
- Equipment code
- Name of the person who cleaned the equipment
- Date of cleaning
- Name of the person who verified the cleaning

Stick the CLEAN label in the appropriate section at the end of this batch record as an evidence that you used a clean equipment

1.4 Labeling

All the containers in contact with the powder need to be properly LABELLED BEFORE the powder is incorporated.
All the printings from printers need to be LABELLED BEFORE the measurement takes place

1.5 SAFETY

Wear gloves and masks at ALL TIMES
BEWARE OF THE BLADE!

1.6 WORKPLACE

YOUR WORKPLACE NEED TO BE CLEANED BEFORE YOU LEAVE THE LABORATORY

Signature prior to the trials Date/sign	Project Manager	QA Manager	Scientists	Signature during the trials Date/sign	Scientist1	Scientist2	Check signature Date/sign	QA Manager
VN 7/04/22	VJ 07 APR 2022.	Yunus		Yunus	K. S			Tanay 07 APR 2022.

7/4/2022 7/4/2022

7/4/2022 7/4/2022

Practical Batch Record- Wet Granulation- V04

Page 3 of 19

DoE Trial No:

CPI

Practical group No:

II

2 PROTOCOL INFORMATION

CAREFULLY read the protocol and this batch record and fill the table below with the appropriate information

2.1 DoE trial

DOE Trial number : CPI(Lower Shear)

Fill the green cells with the appropriate values corresponding to the trial:

factors	Factor 1	Factor 2	Factor 3	Factor 4
Name	Diluent type O	MgSt quantity	Mixer type	NA
Matrix level (-1 or +1)	0	0	-1	V3
Experimental level (real value)	1/2 Lactose 1/2 MCC	1.5%	Low Shear	07 APR 2022.

2.2 Equipment required

Manufacture Equipment	Code No.	Testing equipment	Code No.
1 Top load balance	ARC120	1 Timers	NA. YL
2 Low shear Mixer	213-34	2 Moisture content balance	223-59
3 Food dehydrator	223-51	3 Funnel	223-59
4 Granulator / Screen	66549	4 Tap density tester	13b-2
5 Turbula	060 260	5 Electromagnetic sieve shaker	13b-3
6 Compression machine Shrek	NA	6 Analytical balance	223-100
7		7 Friability tester	13b-12
8	NA TG 74/22	8 Hardness tester	223-58
9		9 Disintegration tester	223-62

Signature prior to the trials Date/sign	Project Manager	QA Manager	Scientists	Signature during the trials Date/sign	Scientist1	Scientist2	Check signature Date/sign	QA Manager
VN 7/4/22	VS 01 APR 2022.	Yunus		Yunus K. J				Yunus 07 APR 2022.

4/7/22 4/7/22

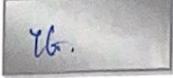
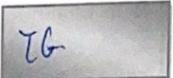
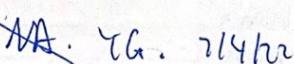
4/7/22 7/4/2022

Practical Batch Record- Wet Granulation- V04

Page 4 of 4

DoE Trial No: CP1

Practical group No: 11

3 DISPENSING		Start Date : 7/4/22 Start Time: 9:45 am
3.1 3.2 3.3	Verify the bubble in the circle (tick the box at the right when OK and put initials) <input checked="" type="checkbox"/>	
	Record the date at which the balance has been calibrated <input checked="" type="checkbox"/>	
	Stick the CLEAN label on the batch record and write with a marker on the plastic bag <ul style="list-style-type: none"> • name of product, <input checked="" type="checkbox"/> • TRIAL batch number <input type="checkbox"/> • date <input type="checkbox"/> • initials <input type="checkbox"/> 	Verified by YG. 7/4/22 9:45 am
	Open the plastic bag, roll the top and put on the top load balance With a dedicated spoon or spatula weigh the excipients and record the amount in the table next page and print the ticket if available. Close the plastic bag and put in the bigger batch plastic bag	
End date : 7/4/22 End time: 10:20 am Done by: VS. 7/4/22	Comments: 	

Signature prior to the trials Date/sign	Project Manager	QA Manager	Scientists	Signature during the trials Date/sign	Scientist1	Scientist2	Check signature Date/sign	QA Manager
VN 7/4/22	VS 07 APR 2022	Yuenho		Yuenho	K. F			07 APR 2022

7/4/22

7/4/22 4/7/2022

7/4/2022

Practical Batch Record- Wet Granulation- V04

Page 5 of 19

DoE Trial No: CP1

Practical group No: 11

Material		Material name		% of the material in the formula (see protocol)	Quantity Required in g	Quantity weighed (recorded by scientist)	Quantity read (recorded by the QA manager)	Checked by (QA manager)
Batch Number	Expiry date	Commercial name	Common name					
10015653	12/22	Pharmatose 200 M	Lactose	2.5	15.00	15.00	15.00	YEG
561021913	06/22	Avicel PH102	Micocrystalline cellulose	32.75	98.25	98.25	98.25	YEG
10015653	12/22	Pharmatose 200M	Lactose	32.15	98.15	98.15	98.15	YEG
41511236W0	04/21	K30	Povidone	4	12.00	11.99	11.99	YEG
M5840	01/23	AC-DI-SOL	Crosmelllose Sodium	4	12.00	12.00	12.00	YEG
18042606	AA. TC 114/22		Magnesium stearate (is weighed after the granule has been made and weighed)	1.5	4.90	4.90	4.90	YEG
			TOTAL	100	300	300.391	300.391	YEG

BEFORE WEIGHING MAKE SURE THAT THE QUANTITIES CORRESPOND TO THE CORRESPONDING TRIAL FORMULA

* WRONG VALUES USED FOR CALCULATION
VS OT APR 2022.

QA signature and date and initials:

Signature prior to the trials Date/sign	Project Manager	QA Manager	Scientists	Scientist1	Scientist2	Check signature Date/sign	QA Manager
				Signature during the trials Date/sign	Signature during the trials Date/sign	Date/sign	Date/sign
VN 7/10/22	VN	VJ 07 APR 2022	YEG 07 APR 2022	YEG 07 APR 2022	YEG 07 APR 2022	YEG 07 APR 2022	YEG 07 APR 2022

7/14/22 7/14/22

Practical Batch Record- Wet Granulation- V04

Page 6 of 8

DoE Trial No:	CP /
Practical group No:	II

MIXER USED : Low Shear

BEWARE THE SHARP BLADE IF YOU USE THE HIGH SHEAR MIXER

4 STEP 1: DRY MIXING Stick CLEAN label <i>V.S.</i> <i>07 APR 2022</i>		Start Date : 7/4/22 Start Time: 10:30 am	
4.1 Put the ingredients in the order described in the mixer bowl IN A SANDWICH FORM, close the lid, put a plastic bag around the mixer to contain the powder and proceed with the mixing:	<ul style="list-style-type: none"> • $\frac{1}{2}$ of $\frac{1}{2}$ or diluent: name of the diluent: <input checked="" type="checkbox"/> • Paracetamol <input checked="" type="checkbox"/> • Binder- name of the binder: <input checked="" type="checkbox"/> • Disintegrant- name of the disintegrant <input checked="" type="checkbox"/> • Remainder diluent: name of diluent: <input checked="" type="checkbox"/> • Speed (only for low shear mixer): <input checked="" type="checkbox"/> • Time: 3 min to be measured with a timer <input checked="" type="checkbox"/> <p>Stop the mixer</p>		
	End date : 7/4/22 End time: 10:45 am Done by: V.S.		
	Comments: <i>NA</i> <i>V.S. 07 APR 2022.</i>		
	Verified by: YH. 7/4/22		

Signature prior to the trials Date/sign	Project Manager	QA Manager	Scientists	Signature during the trials Date/sign	Scientist1	Scientist2	Check signature Date/sign	QA Manager
<i>VN</i> 7/4/22	<i>V.S.</i> 07 APR 2022.	<i>Y.G.</i>		<i>Y.G.</i> 07 APR 2022.	<i>K.S.</i>			<i>Y.H.</i> 07 APR 2022.

7/4/22

7/4/22

7/4/22

Practical Batch Record- Wet Granulation- V04

Page 7 of 19

DoE Trial No: CP 1

Practical group No: 1)

5 STEP 2: GRANULATION		Start Date : 7/4/22 Start Time: 10:45 am
5.1	In a measuring cylinder, measure the correct volume of water. <ul style="list-style-type: none"> Volume of water: <u>40mL</u> <input checked="" type="checkbox"/>	Verified by: <u>YH.</u> 7/4/22
5.2	When ready start the mixer and immediately start to slowly pipette the water into the powder <ul style="list-style-type: none"> Speed (only for low shear mixer): Time of pouring the water : 30s (timer) <input checked="" type="checkbox"/> Continue to mix: <ul style="list-style-type: none"> Time: 2.5 min (total of granulation 3 min) <input checked="" type="checkbox"/> Remove the granule from the mixer and spread it directly in the dryer.	Verified by: <u>YH.</u> 7/4/22
End date : 7/4/22 End time: 10:50 am Done by: <u>VS.</u>		Comments: <i>NA</i> <i>VS</i> <i>07 APR 2022</i>

Signature prior to the trials Date/sign	Project Manager	QA Manager	Scientists	Signature during the trials Date/sign	Scientist1	Scientist2	Check signature Date/sign	QA Manager
VN 7/4/22	VJ 07 APR 2022	YG 07 APR 2022		TG 07 APR 2022	k.r			<i>[Signature]</i> 07 APR 2022

7/4/22

7/4/22

7/4/22

Practical Batch Record- Wet Granulation- V04

Page 8 of

DoE Trial No:

CPI

Practical group No:

II

6 STEP 3: DRYING

Stick CLEAN label

NA-TG

Start Date : 7/4/22

Start Time: 10:55am

6.1	Lay pieces of large mesh cloth on each tray of the dryer and spread the wet granulation the cloth by layer and start the drying <ul style="list-style-type: none"> • Temperature: 75°C <input checked="" type="checkbox"/> • Time: 20 min to start with <input checked="" type="checkbox"/> 		Verified by: YH. 7/4/22
	End date : 7/4/22 End time: 11:15 am Done by: VS.	Comments: First moisture measurement was 3%, therefore dried for 40 mins (20 mins extra).	
6.2	TEST: Moisture content Stick the clean sticker <input checked="" type="checkbox"/> NA-TG. <ul style="list-style-type: none"> • IR balance <input checked="" type="checkbox"/> • Put 2.5g of powder in the tray <input checked="" type="checkbox"/> • Time: 3 min <input checked="" type="checkbox"/> • 110°C <input checked="" type="checkbox"/> • Reading every min <input checked="" type="checkbox"/> 		Result 1.1 %
	End date : 7/4/22 End time: 12:27 pm Done by: VS.	Comments: <i>NA VS 07 APR 2022.</i>	

Signature prior to the trials Date/sign	Project Manager	QA Manager	Scientists	Signature during the trials Date/sign	Scientist1	Scientist2	Check signature Date/sign	QA Manager
VN 7/04/22	VJ 07 APR 2022	YH 07 APR 2022		YH 07 APR 2022	K.J			<i>Luffy</i> 07 APR 2022

7/4/22

7/4/22

7/4/22

Practical Batch Record- Wet Granulation- V04

Page 9 of 19

DoE Trial No: CP1

Practical group No: 11

7 STEP 4: SCREENING

Stick CLEAN label VN YG.

Start Date : 7/4/22

Start Time: 13:26 pm

7.1

Prepare the granulator according to the SOP. Place a large Plastic bag around the bottom of the granulator to recoup the screened granule and screen the granule until no more granule can pass through the mesh

Verified by:

YH.

7/4/22

End date : 7/4/22

End time: 13:31 pm

Done by: VS.

Comments:

NA
VS 07 APR 2022.

8 STEP 5: LUBRICATION

Stick CLEAN label VN YG.

Start Date : 7/4/22

Start Time: 13:34 pm

8.1

TEST: Yield: accurately weigh the dry granule on a top-load balance by taring another plastic bag and putting the bag containing the granules on the balance.

91.7%

Calculate the quantity of lubricant that is needed knowing that the lubricant is 0.5% of the total formula

Equation:

$$\text{Yield} \times 1.5\% \text{ mass of screened granules} \times 1.5\% = \text{mass MgSt required.}$$

- Quantity of Mg Stearate required?

Verified by:

YH.

Verified by:

YH.

5.3202
5.32g

TG . 81412

5.3202
5.32g

TG . 81412
Verified by
YH . 7/4/22

End date : 7/4/22

End time: 13:48 pm

Done by: VS & YG

Comments:

NA
VS 07 APR 2022.

Signature prior to the trials Date/sign	Project Manager	QA Manager	Scientists	Signature during the trials Date/sign	Scientist1	Scientist2	Check signature Date/sign	QA Manager
	VN 7/4/22	VS 07 APR 2022.	YG		TG 7/4/22	K.T 7/4/22		YH 01 APR 2022.

7/4/22

7/4/22

7/4/22

Practical Batch Record- Wet Granulation- V04

Page 10 of

DoE Trial No:

ep1

Practical group No:

II

8.2	Put the dry granule in the glass jar then add the Mg Stearate. Put the glass jar in the Turbula. Start the mixing. • Speed: 72 rpm <input checked="" type="checkbox"/> • Time: 3min <input checked="" type="checkbox"/>		Verified by: YH. 4/4/22 7
	End date : 7/4/22 End time: 14:03 pm Done by: VS.	Comments: NA VS 07 APR 2022	
8.3	TEST: Flowability: weigh 50 g of lubricated powder. Put the funnel in a stand above a 250 ml measuring cylinder that is specific to the tap density tester. Stop the hole of the funnel with your thumb and pour the 50 g of powder in the funnel. Remove your thumb at the exact same time as you activate the timer and stop the timer when the 50 g have gone through the funnel. Make sure not to touch the cylinder if possible when determining the volume. • Flow time in s (no decimal places)	Weight 50.15 g	Result 25
	End date : 7/4/22 End time: 14:40 Done by: VS. & SK.	Comments: NA VS 07 APR 2022	

Signature prior to the trials Date/sign	Project Manager	QA Manager	Scientists	Signature during the trials Date/sign	Scientist1	Scientist2	Check signature Date/sign	QA Manager
VN 7/04/22	VS 07 APR 2022.	TG.		YG. 07 APR 2022.	K.S			Willy 07 APR 2022.

7/4/22

7/4/22

7/9/22

Practical Batch Record- Wet Granulation- V04

Page 11 of 19

DoE Trial No: CPI
Practical group No: 11

8.4	TEST: tapped density: record the bulk volume in the measuring cylinder					Result	
	Cover the measuring cylinder with a metal foil, but not covering the graduated lines on the cylinder. Fit the measuring cylinder in the tapped density tester and set 10 taps according to USP 2. Record the volume after 10 taps and record the powder volume in the table below.	Set the machine for having the following number of taps records and record the volumes in the table below	Pour the powder in the mesh for particle sizing and then clean the cylinder with a cloth (mesh) and pen.				
	Volume (ml)	V0	V10	V500	V1250	(V2500)	
		116	94	90	85	85	
	Density (g/ml)	0.43	0.53	0.56	0.59	0.59	
	Calculate the Carr index:	$C = 100 \frac{V_{10} - V_{500}}{V_{10}}$					
	End date: 4/7/22	Comments: NIA, VN 7/04/22					
	End time: 16:01pm						
	Done by: VS & VN.						
	Sig.						

Signature prior to the trials Date/sign	Project Manager	QA Manager	Scientists	Signature during the trials Date/sign	Scientist1	Scientist2	Check signature Date/sign	QA Manager
VN 7/04/22	VJ 07 APR 2022	TG 7/14/22		TG 7/14/22	K.J			LAWI 07 APR 2022

7/14/22 7/4/22

Practical Batch Record- Wet Granulation- V04

DoE Trial No: CP1
 Practical group No: 11

8.5	<p>TEST: particle size: Measure the weight of each screen:</p> <table border="1" style="margin-left: 20px;"> <tr> <th></th> <th>bottom</th> <th>125µm</th> <th>180 µm</th> <th>355 µm</th> </tr> <tr> <td>Mass in g</td> <td>282.16</td> <td>299.60</td> <td>297.61</td> <td>314.38</td> </tr> </table> <p>Set up the sieve shaker and pour the tapped powder into the top screen and run the test:</p> <ul style="list-style-type: none"> • Time: 10 min <input checked="" type="checkbox"/> • Vibration: 20 <input checked="" type="checkbox"/> <p>Record the weights of each screen containing the powder and KEEP each fraction in separate labeled plastic bags;</p> <table border="1" style="margin-left: 20px;"> <tr> <th></th> <th>bottom</th> <th>125 µm</th> <th>180 µm</th> <th>355 µm</th> </tr> <tr> <td>Mass in g</td> <td>302.10</td> <td>308.90</td> <td>310.61</td> <td>321.12</td> </tr> </table> <p>Calculate the difference and the percentage of powder in each particle size:</p> <table border="1" style="margin-left: 20px;"> <tr> <th></th> <th>bottom</th> <th>125 µm</th> <th>180 µm</th> <th>355 µm</th> </tr> <tr> <td>%</td> <td>40.7</td> <td>19.0</td> <td>26.5</td> <td>13.8</td> </tr> </table> <p>Clean cylinder with mesh and pen before handing it to other group.</p>						bottom	125µm	180 µm	355 µm	Mass in g	282.16	299.60	297.61	314.38		bottom	125 µm	180 µm	355 µm	Mass in g	302.10	308.90	310.61	321.12		bottom	125 µm	180 µm	355 µm	%	40.7	19.0	26.5	13.8	<p>VERIFIED BY <i>[Signature]</i> 08 APR 2022.</p>
	bottom	125µm	180 µm	355 µm																																
Mass in g	282.16	299.60	297.61	314.38																																
	bottom	125 µm	180 µm	355 µm																																
Mass in g	302.10	308.90	310.61	321.12																																
	bottom	125 µm	180 µm	355 µm																																
%	40.7	19.0	26.5	13.8																																
	<p>End date : 08/04/2022 End time: 11:42 a.m. Done by: VS.</p>	<p>Comments: <i>[Signature]</i> VS 08 APR 2022.</p>																																		

9 STEP 6: TABLETTING Stick CLEAN label			Start Date : 7/4/22 Start Time: 3:50 pm				
9.1	<p>Proceed with tabletting on the Shrek machine. Produce at least 200 tablets at the set COMPRESSION FORCE.</p> <p>Stop the machine when the powder volume is low.</p>					Verified by: <i>VS</i> 08 APR 2022.	
	<p>End date : 7/4/22 End time: 4:10pm Done by: SK.</p>	<p>Comments: <i>[Signature]</i> VS 7/4/22</p>					

Signature prior to the trials Date/sign	Project Manager	QA Manager	Scientists	Signature during the trials Date/sign	Scientist1	Scientist2	Check signature Date/sign	QA Manager
	VN	VS	TG.	TG.	K. T			<i>[Signature]</i> 08 APR 2022.

8/4/22

8/4/22

8/4/22

Practical Batch Record- Wet Granulation- V04

Page 13 of 19

DoE Trial No: CP1

Practical group No: 11

9.2	TEST: Mass/Hardness: On an analytical balance, weigh each 10 tablets individually and place them back in the well plate. Measure the hardness and thickness of the tablets in the tablet holder in order to get the mass, the hardness and the thickness of each individual tablet.						Verified by: VS 8/4/22	
	Mass (balance with 4 decimal places) in mg	Tablet 1	Tablet 2	Tablet 3	Tablet 4	Tablet 5		
	0.5118	0.5079	0.5316	0.5032	0.5015			
	Hardness (hardness tester) (N)	24.799	25.449	21.652	22.524	19.242		
	Thickness (Caliper) (mm)	N/A, VN 8/04/22 not required as per video protocol						
	Tablet 6	Tablet 7	Tablet 8	Tablet 9	Tablet 10			
	Mass (balance with 4 decimal places) in mg	0.5165	0.5098	0.5116	0.5076	0.5149		
9.3	Hardness (N)	23.935	20.409	18.922	20.073	15.997		
	Thickness (mm)	not required as per video protocol N/A, VN 8/04/22						
	End date : 8/04/22	Comments:						
	End time: 10:55 am	N/A						
	Done by: VN	VN 8/04/22						
9.3	TEST: Disintegration time: place 3 tablets in the basket holes in 2 baskets and place the disks on top of each tablet. Start the disintegration test according to the pharmacopeia. Record the time at which the tablet is disintegrated (small lumps are considered disintegrated) –						Verified by: VS 8/4/22	
	Tablet 1	Tablet 2	Tablet 3	Tablet 4	Tablet 5	Tablet 6		
	Disintegration time (s)	1:16	1:16	1:16	1:16	1:16		
	End date : 8/04/22	Comments:						
	End time: 10:40 am	N/A						
	Done by: VN	VN 8/04/22						

Signature prior to the trials Date/sign	Project Manager	QA Manager	Scientists	Signature during the trials Date/sign	Scientist1	Scientist2	Check signature Date/sign	QA Manager
8/04	VN	VS	YG.	YG	K.S			harry

8/4/2022

8/4/2022

8/4/2022

8/4/2022

Practical Batch Record- Wet Granulation- V04

Page 14

DoE Trial No:

CPI

Practical group No:

II

9.4	<p>TEST: friability:</p> <ul style="list-style-type: none"> accurately weigh approximately 6.5 g of tablets (4 g if average weight of individual tablet is more than 650 mg) <input checked="" type="checkbox"/> Run the test at 100 rotations (type in 100) <input checked="" type="checkbox"/> <p>Remove the tablets from the wheel, dedust them and weigh them again. Calculate the friability in % = (before-after)/before x 100</p> <table border="1"> <thead> <tr> <th></th><th>Before</th><th>After</th><th>Friability %</th></tr> </thead> <tbody> <tr> <td>Mass (g)</td><td>6.6061</td><td>6.1827</td><td>6.4%</td></tr> </tbody> </table>		Before	After	Friability %	Mass (g)	6.6061	6.1827	6.4%	Verified by: YH. 08/04/2022
	Before	After	Friability %							
Mass (g)	6.6061	6.1827	6.4%							
	<p>End date: 08/04/2022 End time: 11:08 Done by: VS.</p>	<p>Comments: NA VS 08 APR 2022.</p>								

Signature prior to the trials Date/sign	Project Manager	QA Manager	Scientists	Signature during the trials Date/sign	Scientist1	Scientist2	Check signature Date/sign	QA Manager
VN 8/104	V3 08 APR 2022.	YG		YG	K. S			Wifly 08 APR 2022.

8/4/22

8/4/22 5/4/2022

Practical Batch Record- Wet Granulation- V04

Page 15 of 19

DoE Trial No: 1P1

Practical group No: 11

10 STEP 7: TABLETTING GAMLEN

Stick CLEAN label

Start Date :
Start Time:

Verified by:

10.1

On the Gamlen tablet press, produce 4 tablets with the granule

- Mass of each tablet: 200mg
- Compaction force: 500kg

Test: hardness of 4 tablets from the whole granule :

Measure the corresponding hardness and tensile strength. Record the force versus displacement curve and the ejection force and calculate tensile strength

	Tablet 1	Tablet 2	Tablet 3	Tablet 4
Mass (mg)				
Hardness (N)				
Tensile strength (MPa)				
FORCE/ DISPLACEMENT recorded				
EJECTION FORCE				

NA VS 08 APR 2022.

10.2

On the Gamlen tablet press, produce 4 tablets with the fractions of PARTICLE SIZE TEST

- Mass of each tablet: 200mg
- Compaction force: 500kg

Test: hardness of 4 tablets from the whole granule :

Measure the corresponding hardness and tensile strength. Record the force versus displacement curve and the ejection force and calculate tensile strength

	Tablet 1 BOTTOM	Tablet 2 125microns	Tablet 3 180 microns	Tablet 4 355 microns
Mass (mg)				
Hardness (N)				
Tensile strength (MPa)				
FORCE/ DISPLACEMENT recorded				
EJECTION FORCE				

* NOT REQUIRED AS PER VIDEO

Protocol

VJ DJ APR 2022.

Signature prior to the trials Date/sign	Project Manager	QA Manager	Scientists	Signature during the trials Date/sign	Scientist1	Scientist2	Check signature Date/sign	QA Manager
			VS					

DJ APR 2022.

Practical Batch Record- Wet Granulation- V04

Page 16

DoE Trial No:

CPI

Practical group No:

II

End date :	Comments:
End time:	
Done by:	

NA
VS 08 APR 2022.

* NOT REQUIRED AS PER VIDEO PROTOCOL.

11 GENERAL COMMENTS ON THE BATCH

NA
VS 08 APR 2022.

Signature prior to the trials Date/sign	Project Manager	QA Manager	Scientists	Signature during the trials Date/sign	Scientist1	Scientist2	Check signature Date/sign	QA Manager
		VJ		VS 08 APR 2022.				Willy 08 APR 2022.

Practical Batch Record- Wet Granulation- V04

Page 17 of 19

DoE Trial No:

c96

Practical group No:

11

12 CLEAN LABELS

CLEAN LABELS:

NA
VS
08 APR 2022.

Signature prior to the trials Date/sign	Project Manager	QA Manager	Scientists	Signature during the trials Date/sign	Scientist1	Scientist2	Check signature Date/sign	QA Manager
				NA VS 08 APR 2022.				

08 APR 2022.

Practical Batch Record- Wet Granulation- V04

Page 18

DoE Trial No: CP

Practical group No: 11

~~NA~~
YJ
01 APR 2022.

Signature prior to the trials Date/sign	Project Manager	QA Manager	Scientists	Signature during the trials Date/sign	Scientist1	Scientist2	Check signature Date/sign	QA Manager

Practical Batch Record- Wet Granulation- v04

Page 19 of 19

DoE Trial No:

CP

Practical group No:

R1

NA
KJ
0d APR 2022.

Signature prior to the trials Date/sign	Project Manager	QA Manager	Scientists	Signature during the trials Date/sign	Scientist1	Scientist2	Check signature Date/sign	QA Manager